"Thínk Bíochemístry! Thínk Vector..."





UNICORN 230

Random Access Biochemistry Analyzer

200 Tests/Hour, Double Reagent

- Fully Automatic, Discrete, Random access Biochemistry analyzer.
- Throughput: 200 test/hour with STAT function.
- Integrated fully closed optical system.
- 44 Semi Permanent pyrex cuvettes and on board washing system.
- · On board Refrigeration of reagent & sample compartment
- · Available with ready to use system packs.

Multifunctional precision sampling probe

- Polished internally and externally, effective probe washing internally and externally.
- Liquid detection sensor, Accurate detection of reagents and sample volume.
- Collision Protection and automatic adjustment of probe.



Excellent mixer design

- Teflon coated mixer, Prevents water dropping (reducing carry-over)
- Homogenous mixing effect with standard mixing procedure steps.



High precision ceramic syringe

- Durable ceramic piston
- Accurate dispensing with step of 0.1µl
- Effective Pipetting at every dispension.
- Air purge for removing bubble



Reagent & sample tray

- 40 reagent positions (40 for R1 & 40 for R2)
- 40 sample positions; micro cup & test tube can be used
- 24 hours cooling system using peltier pad inside,
- STAT Position possible at any sample position
- 40 sample Positions, QC, Calibrators can be load at any Places.



Reaction tray

- 44 reaction pyrex cuvettes, Automated on board7 steps washing laundry
- Reaction volume as low as 150 μl
- Stable & accurate temperature (37±0.1°C) for reaction
- Separate power switch for refrigeration
- Window for easy changing of lamp

Unicorn 230 System Packs


- Barcoded ready to use liquid stable
- High on board stability
- Total solution of reagents, calibrators and controls

TOTOGRAFIADO. 40 MI

Advanced operating system

- User Friendly Iconic based operation menu, STAT Facility
- Administration Test Profile Panel, Easy Analysis of Data, Option of calculated test programming, Customization of printing formats
- Sample, Reagent and Reaction Disc Status, Sample Blank and Reagent Blank Auto Calculation, QC and Calibration Status, Reagent Volume Monitoring, Substrate Depletion Detection, Abnormal Status Warning, Automated features for Dilution (Pre and Post) and Rerun Option, Real time curve.
- Windows XP, Windows 7, Windows 8, Windows 10
- TCP/IP network
- 6 Optional Formats, Support Customization



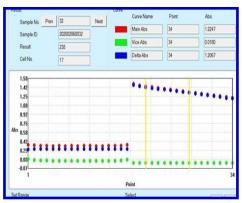

UNICORN 230

Random Access Biochemistry Analyzer

Technical Specifications

General Information:

Machine Type : Random access

Test Speed : Constant 200 tests/hour

(Mono Reagent or Double Reagent)

Test Principle : Colorimetric method, turbidimetry

Immuno turbidimetry

Test Method : 1 Point End, 2 Point End, Fixed Time

and Kinetics

Calibration Type: Linear and Non Linear

Sample Unit:

Sample Tray : 40 Sample Positions

Sample Tube : Micro Cup & Test Tube & Blood Collection

Tube (Φ 12 - 13) mm * (25 - 100) mm

Sample Volume : 2 - 30μL, Step by 0.1 μL

Sample Probe : Liquid surface detection, timely tracking

for liquid volume. Collision Protection

Sample Barcode: Sample barcode reader (Option)

Reagent Unit:

: 40 Reagent Positions (40 for R1 or 40 Reagent Tray

for R2) 24 hours cooling system with

paltier pad inside.

Reagent Volume : $20 - 300 \mu L$, Step by 1 μL

Reaction Unit:

Reaction Tray : 44 reaction new pyrex cuvettes, perfect

transmittance.

Reaction Volume : 150µL - 330µL

Mixing system

Mixer : Paddle-type mixing, effective and

reducing carryover

Optical System:

Light Source : Halogen Lamp

: 340nm; 405nm; 505nm; 546nm; 578nm; Wave Length

630nm; 700nm (4 more options)

Absorbance Range : 0 - 4.0 Abs Resolution : 0.0001 Abs

Operating system

Operating System : Windows XP, Windows 7, Windows 8,

Windows 10

Report Printing : 6 Formats Optional, Support Custom

Formats

Interface : Pattern, TCP/IP network Interface

Working Conditions

Power Supply : AC 110/220V ± 10%, 50/60 Hz, 1000VA

Water consumption: 8L/hour (maximum)

Dimension : 750mm * 550nm * 490mm (L*W*H)

Weight : 70 kgs

Power Supply

110/220V ± 10%, 50/60 Hz, 1000VA, Earthing voltage < 5V online ups (230V) with isolated transformer 2KVA





Exclusively Marketed in India by:

VECTOR BIOTEK PVT. LTD.

Factory & Registered Office: 425, NEW GIDC, KABILPORE, NAVSARI - 396 424. (GUJ.) INDIA Toll free:1800-270-1027 I Tel:+91-93757 94103 I Fax:+91-2637-265105 E-mail: service@beaconindia.com I Web: www.vectorbiotekindia.com CIN No. U33112GJ2002PTC040978

/BPL/JULY20/UNI230/Ver1.0/1000